

ViroQ SARS-FluA/B-RSV RT-PCR Kit

COMING in October

Your benefits

One-step multiplex real time SARS-CoV-2, Influenza and RSV Test

The ViroQ SARS-FluA/B-RSV - Kit is designed for specific detection of viral RNA genes N and RdRP of SARS-CoV-2, Influenza A and separate B as well as RSV in respiratory samples.

Be confident of negative results

- False negative results due to a lack of cell material in the sample are excluded by simultaneous detection of a universally expressed human gene
- A mix of human RNA and SARS-CoV-2, Influenza and RSV reference RNA is used as positive control (included in the kit)

Our product features - Your benefits

- Highly specific and exclusive detection of the SARS-CoV-2 RdRP and N gene,
 Influenza A, Influenza B and respiratory syncytial virus
- One-Step real time PCR assay easy to set up
- Results in approximately 60 minutes save time
- Works with common RNA-Isolation kits be independent from single suppliers
- Ambient shipping condition no extra costs and handling of dry ice
- Quality made in Germany high supply chain security

Compatible with most real time cyclers

- The multiplex real time PCR kit uses five colours. Four for the specific genes and a cycler-specific colour for the detection of the universally expressed human gene, so that the test is compatible with most common real time cyclers on the market.
- Validated on QuantStudio 6 (Thermo Fischer Scientific), CFX 96[™] (Bio-Rad)
 Lightcycler® 480 (Roche)




Workflow and Interpretation



Example for positive result

RdRP and N gene; Flu A; Flu B; RSV Included Separate Internal control for: CFX / LightCycler® 480 / QuantStudio™

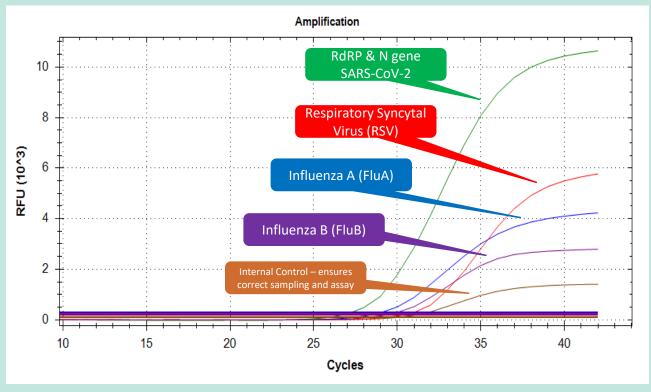

Artificial sample to show all Signals; this case has an extremely low probability of occurrence



Ordering Information

Combination test - for SARS-CoV-2, FluA, FluB, RSV

Cat. No.	Kit	Size	CING
728267	ViroQ SARS-FluA/B-RSV	96 tests	COMING in
728268	ViroQ SARS-FluA/B-RSV	480 tests	October

SARS-CoV-2 test

Cat. No.	Kit	Size
728250	ViroQ SARS-CoV-2	96 tests
728251	ViroQ SARS-CoV-2	480 tests
728263	ViroQ Rapid SARS-CoV-2	96 tests
728264	ViroQ Rapid SARS-CoV-2	480 tests

Ordering:

Phone: +49 (0) 6404 925 450 Fax: +49 (0) 6404 925 460

Email: order@bag-diagnostics.com

Customer Service:

Phone: +49 (0) 6404 925 125 Fax: +49 (0) 6404 925 421

Email: info@bag-diagnostics.com